

Recreating Young's Double Slit Experiment

Ed Lipchus*

Department of Natural Science, Hampshire College, Amherst, MA 01002

Danika Luntz-Martin[†]

Department of Physics, Smith College, Northampton, MA 01063

(Dated: December 9, 2014)

Abstract

I. INTRODUCTION

The double slit experiment performed by Thomas Young was seminal because it demonstrated the wave nature of light.

II. METHODS

We used the TeachSpin Two-Slit Interference, One Photon at a Time (TWS1-B) which is a apparatus designed to perform Young's double slit experiment. The apparatus has two light sources, a 670 nm diode laser and a light bulb with a removable green light filter. There are four slit holders spaced throughout the length of the apparatus, see Figure 1. The first slit, the source slit, is eliminates excess light from the light source and approximates a plane wave. Next is the double slit, our double slit had 0.1mm wide slits and center to center separation of 0.353mm. Following the double slit was a wide single slit attached to a micrometer, called the slit blocker. The slit blocker could be moved to block either or both slits of the double slit. Finally there was a detector slit which allowed us to measure the intensity of light (using the laser) or the number of photons (using the light bulb) over a small region of the x-axis. The micrometer attached to the detector slit allowed us to scan along the x-axis. When we worked with the laser as our light source we used the built in photodetector which gave an output voltage that was proportionate to the light intensity. When we used light bulb as our light source we used the built in photo-multiplier tube (PMT) attached to the TeachSpin Pulse Counter / Interval Timer (PCIT1) to count the number of photons.

We began by aligning the slits within our apparatus. We started by using the laser roughly align the source slit and double slit and ensure that the double slit was not tilted. Then we made fine adjustments to the alignment using the light bulb. Throughout our experiments we used a midrange light bulb power setting of 5. We used a green filter on our light bulb to lower the intensity of the light and narrow the range of wavelengths to 541 - 551nm.

Next we optimized the signal-to-noise ratio for our PMT. To do so we held the PMT voltage constant and varied the discriminator threshold. We found that dial value of the PMT voltage was not consistent with the measured value of the PMT voltage using a digital

Two Slit Apparatus

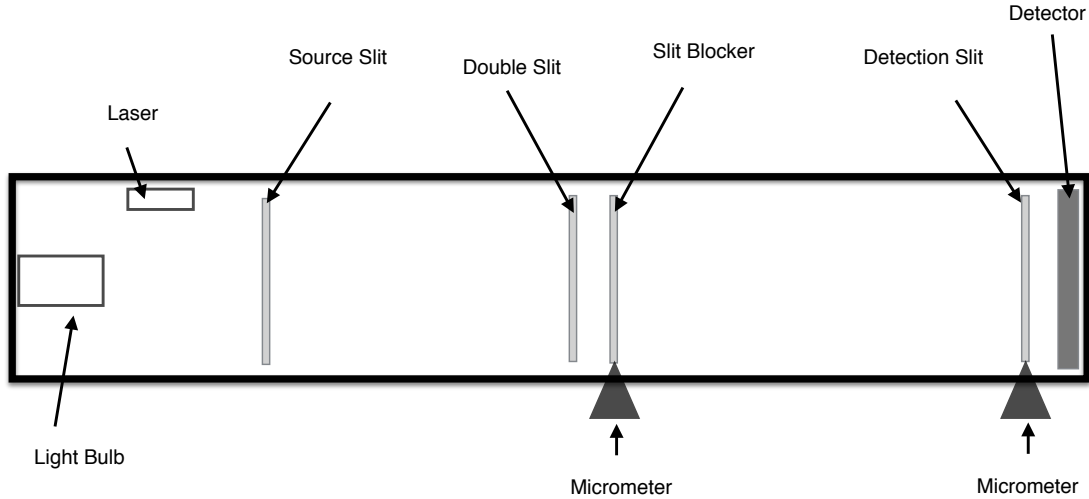


FIG. 1. The experimental set-up of the apparatus we used. The laser can slide into place in front of the light bulb. The source slit is a wide slit that blocks extraneous light from the light source. The slit blocker is another wide slit that can be moved, using the micrometer, to block one of the double slits, effectively turning it into a single slit. The detector slit allowed us to count the number of photon reaching the detector over a small region. The micrometer allowed us to scan across the detector and get a photon count at each x value.

multimeter, so we used the multimeter reading to determine our voltage. For each value of the discriminator threshold we measured four values for our signal plus our noise and four values of just noise with no signal. We then averaged each set of four values and took the ratio to calculate the signal-to-noise. In our case the ratio we actually calculated was signal and noise to noise, this difference should scale our results but not change the overall outcome. The calculated signal-to-noise ratios can be seen for various PMT voltages in Figure ??.

The highest values for the PMT voltage, 900 V and 800 V, had a significantly lower signal to noise ratio. The lower values for the PMT voltage, 600 V and 650 V, required small

values for the discriminator threshold and did not have as high a signal-to-noise ratio. The highest signal-to-noise ratio we obtained was for PMT voltage of 750 V and a discriminator threshold of 15.

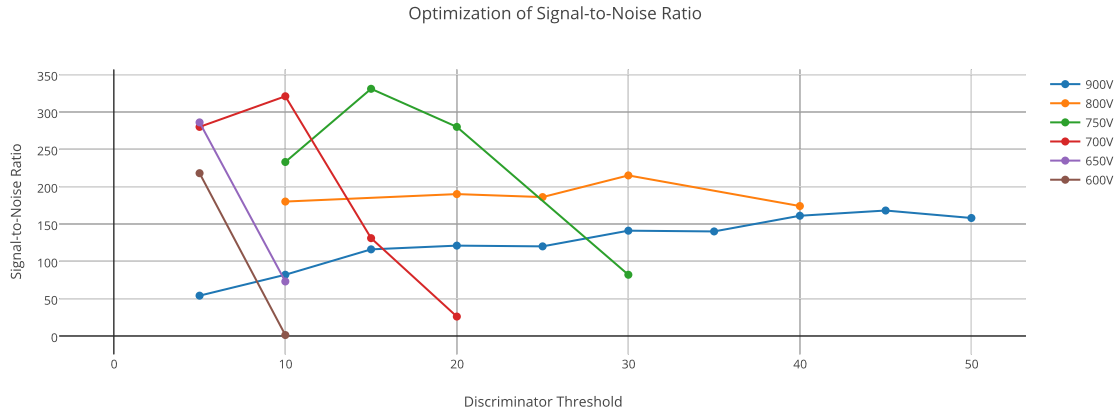


FIG. 2. The signal-to-noise ratios we obtained for a range of discriminator thresholds for various PMT voltages. For the high PMT voltages, 900 V (blue) and 800 V (orange) the signal-to-noise ratio was much lower. For lower PMT voltages, 650 V (purple) and 600 V (brown), the signal-to-noise ratio was again lower. The optimal signal-to-noise ratio was at a PMT voltage of 750 V (green) and a discriminator threshold value of 15.

To collect data we used a Labview program called “TeachSpin Counter V3” that connected to the pulse counter

III. RESULTS

IV. ANALYSIS

V. DISCUSSION

VI. CONCLUSION

* ejl13@hampshire.edu

[†] dluntzma@smith.edu

¹ Jonathan F. Reichert, *TeachSpin Instruction Manuals: Two-Slit Interference, One Photon at a Time (TWS1 - B), Pulse Counter / Interval Timer (PCIT1)* Rev. 1.0, (2013)